



PHOTO TIP SHEET

Send photos we can help!

Where do I send my photos?

- Send all digital photos to quoterequest@cjtmounting.com

IMPORTANT: DO NOT take photos from a standing position. The most helpful photos are taken from where your eyes are level with the seat of the chair. (See Figure 1 on page 2).

What is the best way to take the photos?

- Before taking the photos remove all items from the chair that could block us from seeing a clear view of the chair, bed or walker. (i.e. seatbelts, backpacks, blankets, cup holders etc).
- If possible, place the chair against a solid, light colored back ground.
- Take the photos in a well lighted area or use a flash.
- Cell phone photos are not typically accepted due to the poor photo quality.

What type of photos should I take?

1. Take front and side photos of the person in the chair while another person is holding the device in the desired position. (If you do not have the device use a book or other object that can represent the device).
2. Take photos of each side of the chair.
3. Take photos of the front and the back of the chair.
4. If the chair tilts take photos of the chair in a tilted and upright position.
5. Photos should show a clear view of the entire chair from top to bottom.
6. Take a few closer photos of the seating system tubing and other areas you think may be possible mounting locations. Please, remember not too close or the perspective will be lost.
7. If the client uses a lap tray please be sure to take photos with the lap tray on.
8. Take photos of the brake in the open position and the locked position.

C.J.T. Enterprises, Inc.

P.O. Box 10028
Costa Mesa, Ca 92627




Phone: 714-751-6295
Fax: 714-751-5775

www.cjtmounting.com

quoterequest@cjtmounting.com



PHOTO TIP SHEET

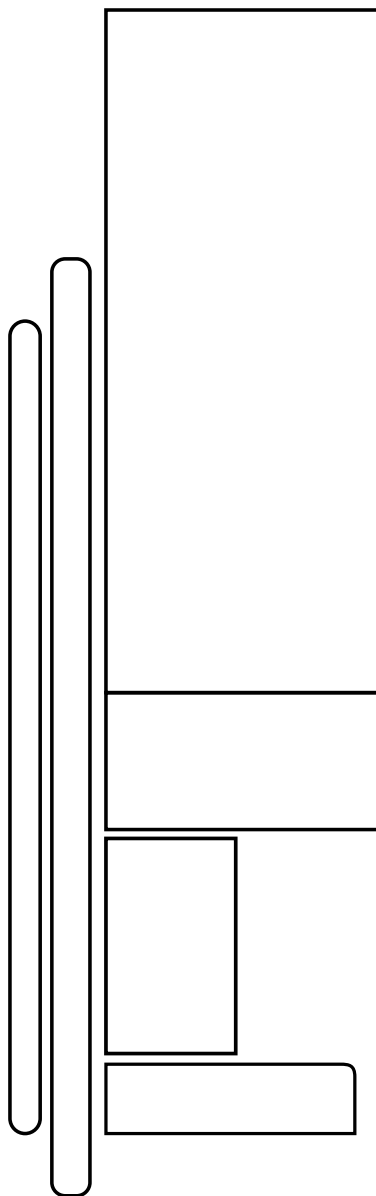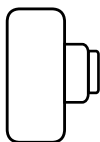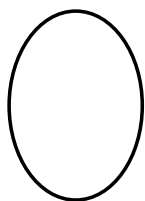

C.J.T. Enterprises, Inc.

P.O. Box 10028
Costa Mesa, Ca 92627

Phone: 714-751-6295
Fax: 714-751-5775

www.cjtmounting.com